


INSTRUCTIONS FOR USE
AIMOVIG® [AIM-oh-vig]
(erenumab-aooe)
injection, for subcutaneous use
70 mg/mL
single-dose prefilled SureClick® autoinjector

This Instructions for Use contains information on how to inject AIMOVIG with a SureClick autoinjector.

If your healthcare provider decides that you or a caregiver may be able to give your injections of AIMOVIG at home, you should receive training on the right way to prepare and inject AIMOVIG. Do not try to inject yourself until you have been shown the right way to give the injections by your healthcare provider or nurse.

The medicine in the AIMOVIG autoinjector is for injection under the skin (subcutaneous injection). See the AIMOVIG Patient Information for information about AIMOVIG.

Getting to know your prefilled autoinjector

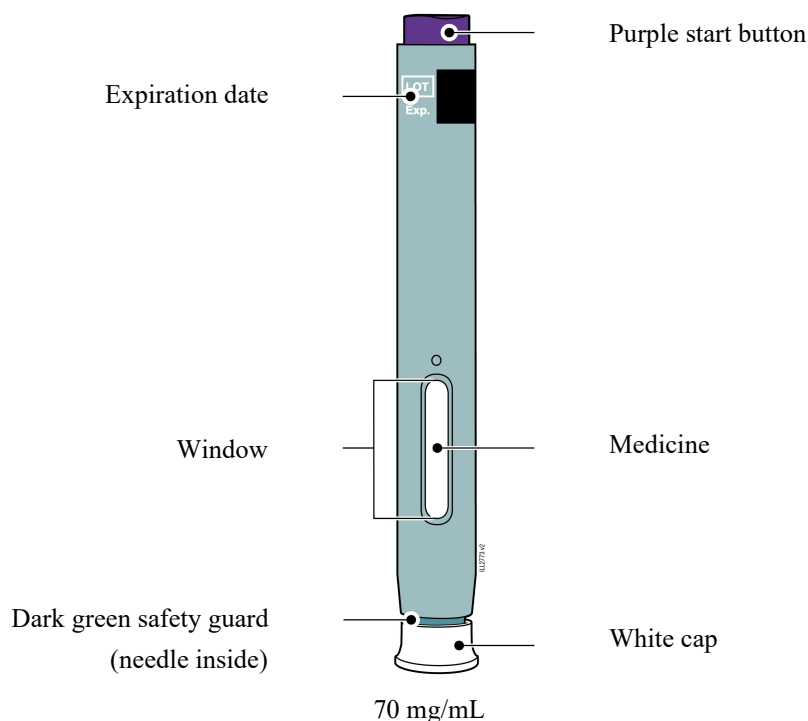

Important Information You Need to Know Before Injecting AIMOVIG

- It is important that you do not try to give the injection until you have fully read and understood this Instructions for Use.
 - Check the autoinjector label and prescription to make sure you have the correct medicine and dose.
 - **Do not** use the autoinjector if the carton is damaged or the seal is broken.
 - **Do not** use the autoinjector after the expiration date on the label.
 - **Do not** shake the autoinjector.
 - **Do not** remove the white cap from the autoinjector until you are ready to inject.
 - **Do not** use the autoinjector if it has been frozen.
 - **Do not** use the autoinjector if it has been dropped on a hard surface. Part of the autoinjector may be broken even if you cannot see the break. Use a new autoinjector, and call 1-800-77-AMGEN (1-800-772-6436).
 - The autoinjector is not made with natural rubber latex.
-

Frequently asked questions:


For additional information and answers to frequently asked questions, visit www.aimovig.com.

Where to get help:

If you want more information or help using AIMOVIG:

- Contact your healthcare provider,
- Visit www.aimovig.com, or
- Call 1-800-77-AMGEN (1-800-772-6436)

Storing and Preparing to Inject AIMOVIG

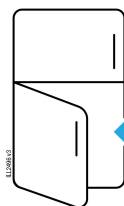


Refrigerate in carton
until ready to use

1 Refrigerate the autoinjector carton until you are ready to use it.

- Keep the autoinjector in the refrigerator between 36°F to 46°F (2°C to 8°C).
- Keep the autoinjector in the original carton to protect it from light or physical damage.
- **Do not** freeze the autoinjector.
- **Do not** store the autoinjector in extreme heat or cold. For example, avoid storing in your vehicle's glove box or trunk.

Important: Keep the autoinjector and all medicines out of the sight and reach of children.

WAIT

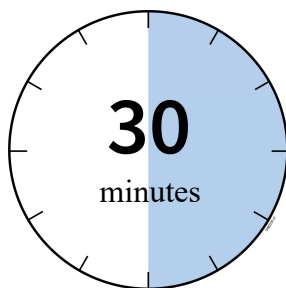

2 Wait 30 minutes for the autoinjector to reach room temperature.


- Let the autoinjector warm up naturally.
 - **Do not** heat the autoinjector with hot water, a microwave, or direct sunlight.
 - **Do not** shake the autoinjector at any time.
 - Using the autoinjector at room temperature makes sure the full dose is delivered and allows for a more comfortable injection.
-



3 You may keep AIMOVIG at room temperature for up to 7 days, if needed.

- For example, when you are traveling, you may keep AIMOVIG at room temperature.
 - Keep it at room temperature between 68°F to 77°F (20°C to 25°C).
 - **Do not** put it back in the refrigerator.
 - Record the date you removed it from the refrigerator and use it within 7 days.


Important: Place the autoinjector in a sharps disposal container if it has reached room temperature and has not been used within 7 days.

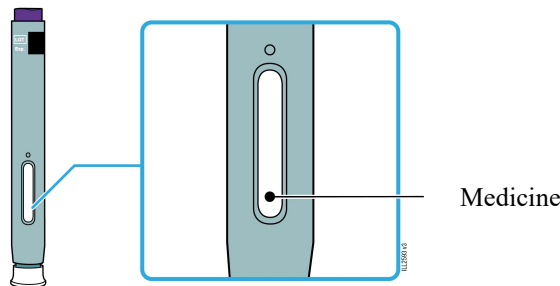

4 Inspect the medicine. It should be clear and colorless to slightly yellow.

- It is okay to see air bubbles.
- **Do not** use the autoinjector if the medicine is cloudy, discolored, or contains flakes or particles.


Important: If the medicine is cloudy, discolored, contains flakes or particles, or if the autoinjector is damaged or expired, call 1-800-77-AMGEN (1-800-772-6436).

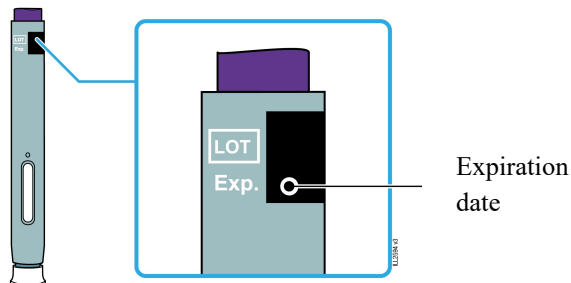

5 Check the expiration date (Exp.) and inspect the autoinjector for damage.


- **Do not** use the autoinjector if the expiration date has passed.
- **Do not** use the autoinjector if:
 - the white cap is missing or loose in carton.
 - it has cracks or broken parts, or
 - it has been dropped on a hard surface.

Getting Ready to Inject AIMOVIG

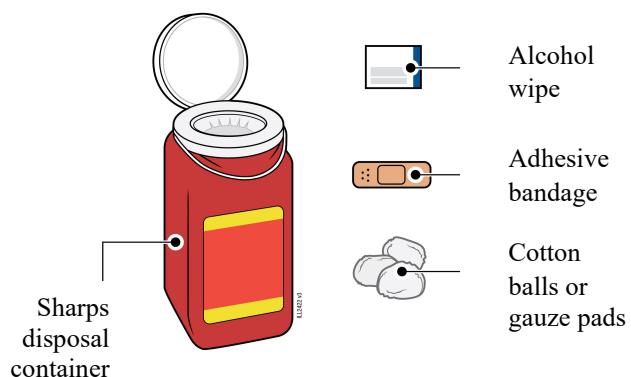

6 Gather and place the following items for your injection on a clean, flat, and well-lit surface:


- AIMOVIG autoinjector (room temperature),
- Sharps disposal container [see Completing the Injection and Disposal],
- Alcohol wipe,
- Adhesive bandage, and
- Cotton balls or gauze pads.

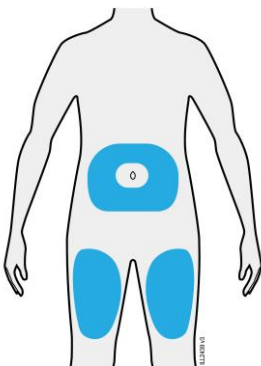

7 Select 1 of these injection sites.

- Select the thigh or stomach (except for 2 inches around the belly button).
- Someone else can inject in your thigh, stomach, or back of the upper arm.

Important: Avoid areas with scars or stretch marks, or where the skin is tender, bruised, red, hard, raised, thick or scaly skin patch, or lesion.

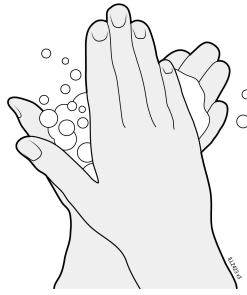

8 Wash your hands thoroughly with soap and water.




9 Clean the injection site with an alcohol wipe.

- Let the skin dry on its own.
 - **Do not** touch this area again before injecting.
-

Injecting AIMOVIG


Important: Only remove the white cap when you can inject right away (within 5 minutes) because the medicine can dry out. **Do not** recap.

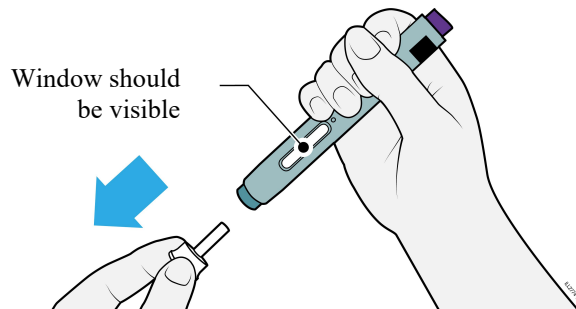

10 Grasp the autoinjector so you can see the window. Pull the white cap straight off. You may need to pull hard.


- **Do not** twist, bend, or wiggle the white cap to pull it off.
 - **Never** put the white cap back on. It may damage the needle.
 - **Do not** put your finger inside the dark green safety guard.
 - It is normal to see a drop of medicine come out of the needle or dark green safety guard.
-

STRETCH

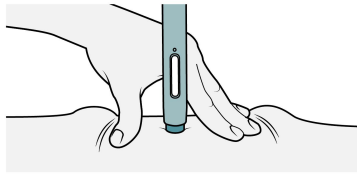

OR


PINCH

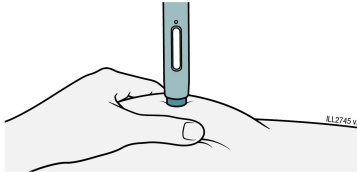

11 Stretch or pinch the skin to create a firm surface at the injection site until the injection is finished. Place the dark green safety guard straight against the skin.

- Make sure you can see the window.
- Make sure the autoinjector is positioned straight on the injection site (at a 90-degree angle).


PUSH
and hold against skin

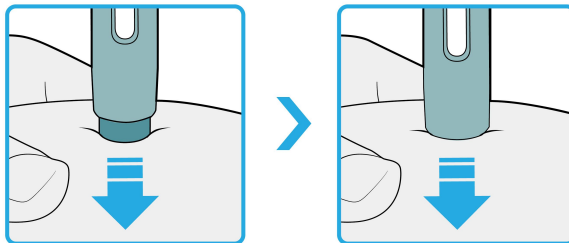

12 Firmly push the autoinjector down until the dark green safety guard stops moving. Hold the autoinjector down, do not lift.

- The dark green safety guard pushes in and unlocks the purple start button.


PRESS
purple start button

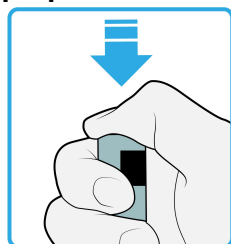

13 Keep pushing the autoinjector down and press the purple start button to start the injection.

- You may hear or feel a click.
- The window starts to turn yellow.
- It is okay to let go of the purple start button.

WATCH
window will turn fully yellow

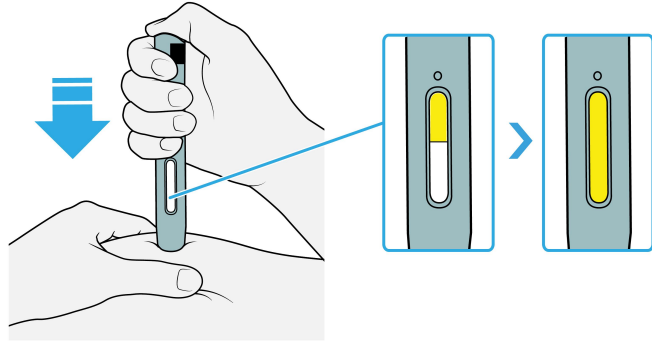


14 Keep pushing the autoinjector down. Wait for the window to turn fully yellow.

- The injection may take up to **15** seconds to complete.
- You may hear or feel a click.
- After the window turns fully yellow, lift the autoinjector away from the skin.
- The dark green safety guard locks around the needle.

Completing the Injection and Disposal

CONFIRM


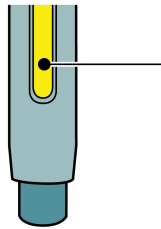


Window is fully yellow

15 Confirm a full dose of medicine was injected.

- **Do not** touch the dark green safety guard.
- A small drop of liquid on the injection site is okay.


Important: If the window has not turned fully yellow or if it looks like the medicine is still coming out, a full dose was not injected. Call your healthcare provider right away.

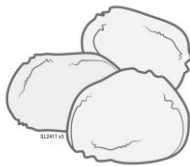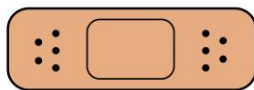

16 Check the injection site.

- **Do not** rub the injection site.
- If there is blood, press a cotton ball or gauze pad on your injection site.
- Apply an adhesive bandage if necessary.

For more information or help call 1-800-77-AMGEN (1-800-772-6436).

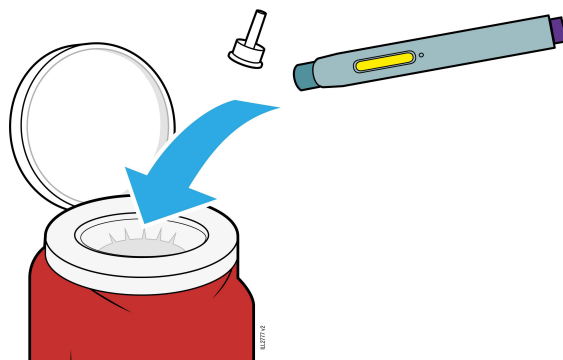

17 Place the used autoinjector and white cap in an FDA-cleared sharps disposal.

Important: Do not throw away the autoinjector in your household trash.

- **Do not** reuse the autoinjector.
- **Do not** touch the dark green safety guard.

Additional information about your sharps disposal container

If you do not have an FDA-cleared sharps disposal container, you may use a household container that is:

- made of a heavy-duty plastic,
- can be closed with a tight-fitting, puncture-resistant lid, without sharps being able to come out,
- upright and stable during use,
- leak-resistant, and
- properly labeled to warn of hazardous waste inside the container.

Disposing of sharps disposal containers:

When your sharps disposal container is almost full, you will need to follow your community guidelines for the right way to dispose of your sharps disposal container.

There may be state or local laws about how you should throw away used needles and syringes.

For more information about safe sharps disposal, and for specific information about sharps disposal in the state that you live in, go to the FDA's website at:

<http://www.fda.gov/safesharpsdisposal>

Do not dispose of your used sharps disposal container in your household trash unless your community guidelines permit this.

Do not recycle your used sharps disposal container.

Keep the autoinjector and sharps disposal container out of the sight and reach of children.

AIMOVIG® (AIM-oh-vig)

AMGEN

Manufactured by:

Amgen Inc.

One Amgen Center Drive

Thousand Oaks, California 91320-1799

U.S. License Number 1080

© 2018 – 2025 Amgen Inc.

All rights reserved.

1XXXXXX

This Instructions for Use has been approved by the U.S. Food and Drug Administration.

Revised: 6/2025 v12